

Did you know that **disorganised**, overloaded baskets can make some instruments **impossible** to clean?



Disorganised baskets can make some instruments impossible to clean because:

- Overloading baskets impede cleaning due to effects like shadowing.
- Unorganised baskets, one layer of instruments on top of the other, make them impossible to clean.
- Flat-lying hollow instruments are prone to biofilm growth because the water stays inside.

But organising baskets takes so much time, and they get scrambled in the machine anyway. Sounds familiar? It's a false economy to think that having a couple more instruments in one load saves time when some of them need to be reprocessed again.

Introducing SmartTray

- Load instruments faster and limit the mess
- Improve cleaning results and limit shadowing
- Extend instrument longevity and limit instrument-on-instrument scratching


Thanks to SmartTray, you load instruments quicker, baskets are typically **lighter**, so **easier** to load into your washer, and instruments are laid out in an **optimal** orientation for cleaning.

Clever solution to your surgical instrument organising troubles!

Options include SmartTray, ScizzoRack, SmartTray MaxFlow Manifolds, and SmartTray LumenRack:

- made from soft polymer and silicone
- withstands temp up to 140°C
- keeps your instruments in place

It's all thought out to make your daily routine more effective, efficient and user-friendly. Order SmartTray today and see the difference for yourself.

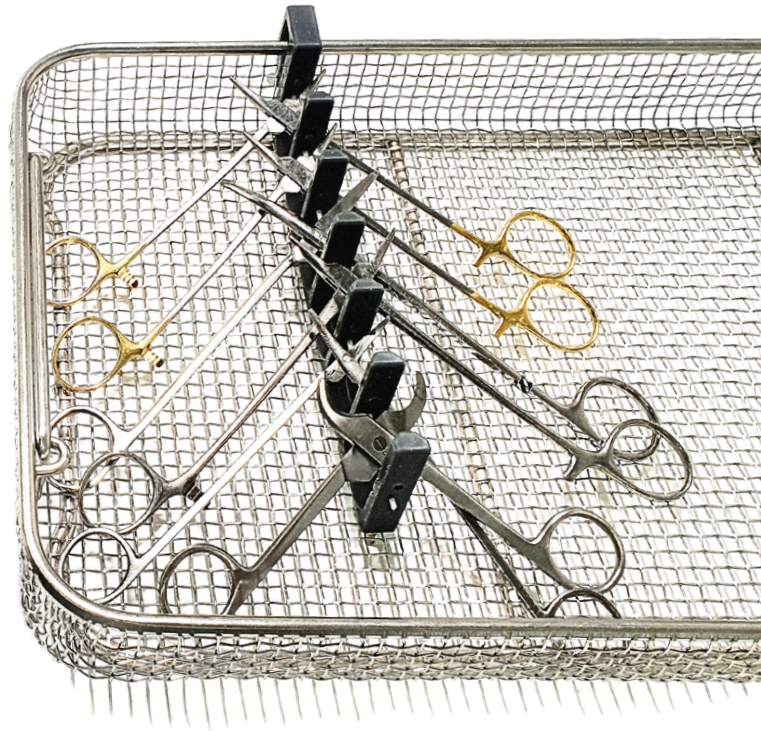

Organise the right way, with SmartTray.

The design ensures your instruments stay in position, and the soft support line gently keeps them in place throughout the cleaning process. What is more, SmartTray's modular design fits instruments of different lengths.


Also, instruments always lay at an angle reducing standing water to the minimum - let gravity work for you! Have your hollow instruments naturally draining so that drying is more effective as water does not reside in the channels, thus limiting the potential for biofilm growth.



ScizzoRack (for 6 and 8 instruments)

Keeps instruments at an angle, making rinsing and drying more effective: it's easier for water to drain, preventing biofilm formation.


Dedicated clip keeps instruments in place during the washing cycle.

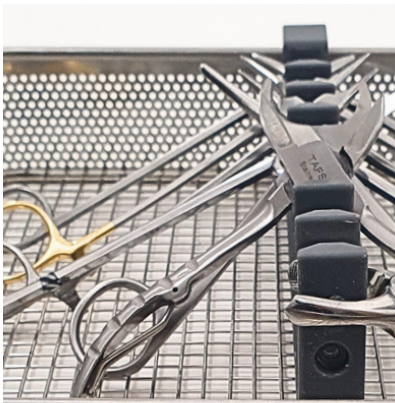


Keeps instruments at an angle, which means no residual water, and faster drying process.



Keeps instruments open, which makes cleaning more effective.

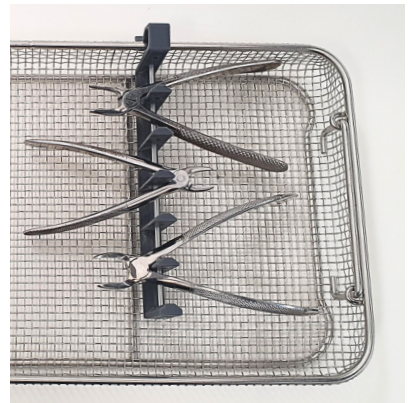

Keeps instruments in their dedicated slots during the washing cycle.

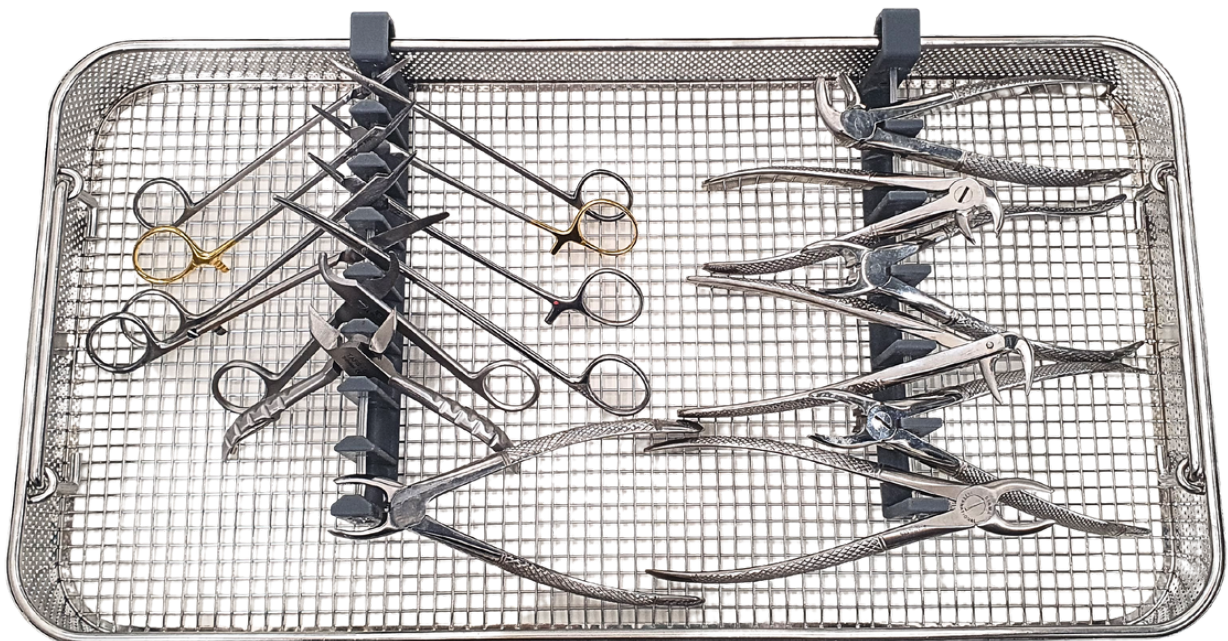

You choose how you want to organise your tray, the design is modular, so clean what needs to be cleaned!

Range SmartTray

Product name	Part no
ScissoRack6 (set of 2)	057.1023
LumenRack6 (set of 2)	057.1024
ScissoRack8 (set of 2)	057.1025
LumenRack8 (set of 2)	057.1026
Manifold MaxFlow6	057.1030
Manifold MaxFlow8	057.1031
Dental 6-instrument Cassette	057.1100
Dental 8-instrument Cassette	057.1101

Go to aseptium.com to find your nearest distributor or email info@aseptium.com for details.